

## SUDHL1 FFPE Cell Pellet

### GENERAL INFORMATION

Product Name:	SUDHL1 FFPE Cell Pellet
Reference Number:	3080-0110 Block 3080-0120 Slide (5µm) 3080-0130 FFPE scroll (20µm)
Date of Manufacturing:	See product label
Lot Number:	See product label
Intended Use:	For research use only

### DESCRIPTION

Cell Line:	SUDHL1
Tissue of Origin:	Lymphoma
Culturing Condition:	RPMI-1640 supplemented with 10% FBS at 37°C with 5% CO <sub>2</sub>
Fixation Condition:	10% neutral buffered formalin (NBF) for 24 hours at 24-27°C
Product Format:	
Block:	Paraffin embedded block. Pellet thickness: ~2mm
Slide:	One unstained section mounted on Superfrost™ Plus slide. Section thickness: 5µm
FFPE Scroll:	One FFPE section in DNase/RNase free tube. Section thickness: 20µm

### SCHEMATICS

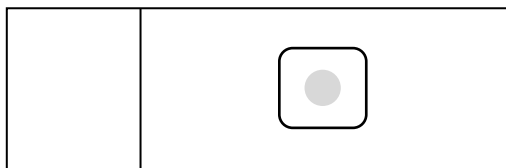


Illustration of an FFPE slide

\*Guideline for general applications, such as immunohistochemistry (IHC) or DNA in situ hybridization (ISH). Certain biomolecules may be less stable during storage

### SAFETY AND PRECAUTIONS

This product does not contain hazardous material. Wear appropriate personal protective equipment (PPE) when handling reagents and biological specimens.

### RECOMMENDED PROCEDURES

#### Staining using FFPE slides:

1. Bake slides at 60°C for 30-60 min.
2. Deparaffinize two times in Xylene or Xylene substitute for 5 min each time.
3. Rinse two times in 100% ethanol for 1 min each time.
4. Air dry slides or rehydrate in ethanol series (95% 2 min, 70% 2 min, 50% 2 min, 1XPBS 2 min).
5. Proceed to staining protocol.

#### Biomolecule extraction from FFPE scrolls:

1. Add 1ml Xylene or Xylene substitute to each tube containing FFPE scrolls and vortex for 10 sec.
2. Centrifuge at full speed for 2 min at room temperature. Remove supernatant without disturbing the pellet.
3. Repeat step 1 and 2.
4. Add 1ml 100% ethanol and mix by vortexing.
5. Centrifuge at full speed for 2 min at room temperature. Remove supernatant without disturbing the pellet.
6. Repeat step 4 and 5.
7. Carefully remove any residual ethanol in the tube without disturbing the pellet.
8. Open the tube and dry at room temperature or up to 37°C for 10min, or until the ethanol has completely evaporated.
9. Proceed to extraction protocol.

### STORAGE AND STABILITY

Storage Condition	2-8°C with desiccation (Recommended)
-------------------	--------------------------------------

### STABILITY

Block	5 years
Scroll	1 year
Slide	1 year

AMSBIO | [www.amsbio.com](http://www.amsbio.com) | [info@amsbio.com](mailto:info@amsbio.com)