

OVCAR3 A-FLX™ FFPE Cell Pellet

GENERAL INFORMATION

| | |
|------------------------|--|
| Product Name: | OVCAR3 A-FLX™ FFPE Cell Pellet |
| Reference Number: | 3090-0710 Block 3090-0720 Slide (5µm) 3090-0730 FFPE scroll (20µm) |
| Date of Manufacturing: | See product label |
| Lot Number: | See product label |
| Intended Use: | For research use only |

DESCRIPTION

| | |
|----------------------|---|
| Cell Line: | OVCAR3 |
| Tissue of Origin: | Ovary |
| Culturing Condition: | DMEM supplemented with 10% FBS at 37°C with 5% CO ₂ |
| Fixation Condition: | 10% neutral buffered formalin (NBF) for 24 hours at 24-27°C |
| Product Format: | |
| Block: | Paraffin embedded block. Pellet thickness: ~2mm |
| Slide: | One unstained section mounted on Superfrost™ Plus slide. Section thickness: 5µm |
| FFPE Scroll: | One FFPE section in DNase/RNase free tube. Section thickness: 20µm |

SCHEMATICS

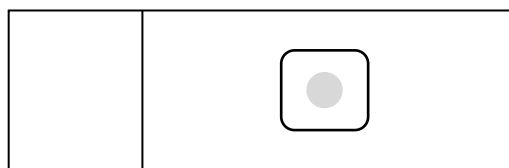


Illustration of an FFPE slide

SAFETY AND PRECAUTIONS

This product does not contain hazardous material. Wear appropriate personal protective equipment (PPE) when handling reagents and biological specimens.

RECOMMENDED PROCEDURES

Staining using FFPE slides:

1. Bake slides at 60°C for 30-60 min.
2. Deparaffinize two times in Xylene or Xylene substitute for 5 min each time.
3. Rinse two times in 100% ethanol for 1 min each time.
4. Air dry slides or rehydrate in ethanol series (95% 2 min, 70% 2 min, 50% 2 min, 1XPBS 2 min).
5. Proceed to staining protocol.

Biomolecule extraction from FFPE scrolls:

1. Add 1ml Xylene or Xylene substitute to each tube containing FFPE scrolls and vortex for 10 sec.
2. Centrifuge at full speed for 2 min at room temperature. Remove supernatant without disturbing the pellet.
3. Repeat step 1 and 2.
4. Add 1ml 100% ethanol and mix by vortexing.
5. Centrifuge at full speed for 2 min at room temperature. Remove supernatant without disturbing the pellet.
6. Repeat step 4 and 5.
7. Carefully remove any residual ethanol in the tube without disturbing the pellet.
8. Open the tube and dry at room temperature or up to 37°C for 10min, or until the ethanol has completely evaporated.
9. Proceed to extraction protocol.

STORAGE AND STABILITY

| Storage Condition | Stability* |
|---------------------------------------|------------|
| Ambient temperature | 6 months |
| 2-8 °C with desiccation (Recommended) | 2 years |
| -20 °C to -10 °C | 5 years |

*Guideline for general applications, such as immunohistochemistry (IHC) or DNA in situ hybridization (ISH). Certain biomolecules may be less stable during storage.

AMSBIO | www.amsbio.com | info@amsbio.com

 **UK & Rest of the World**
184 Park Drive, Milton Park
Abingdon OX14 4SE, UK
T: +44 (0)1235 828 200
F: +44 (0) 1235 820 482

 **North America**
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218

 **Germany**
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880

 **Switzerland**
Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
T: +41(0) 91 604 55 22
F: +41(0) 91 605 17 85