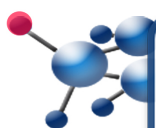


Incomplete Freund's Adjuvant

Catalog # 7002

For Research Use Only - Not Human or Therapeutic Use

DESCRIPTION:	Mineral oil
APPLICATION:	Use as an adjuvant for the immunization of antigen in experimental animals. This incomplete Freund's adjuvant (IFA) is suitable for collagen-induced arthritis (CIA) in rats. This IFA can also be used for a booster injection of type II collagen to increase the incidence of CIA in mice. Note: IFA instead of complete Freund's adjuvant (CFA) should be used to emulsify type II collagen for immunizing rats, since rats are susceptible to adjuvant induced arthritis. For inducing arthritis in rats more efficiently, muramyl dipeptide (MDP) can be added to IFA.
QUANTITY:	5 mL
FORM:	Solution
EMULSION:	Emulsify with an equal volume of collagen using an electric homogenizer. Syringe-syringe method is not recommended to make collagen-IFA emulsion for inducing arthritis, but can be used just for obtaining anti-sera.
STORAGE TEMPERATURE:	4°C in the dark
STABILITY:	1 year
REFERENCES:	William C, et al. J. Exp. Med. 162: 1878-1891 (1985) Trentham D, et al. J. Exp. Med. 146: 857-868 (1977) Taurog JD, et al. J. Exp. Med. 162: 962-978 (1985) Koga T, et al. Arth. Rheum. 23: 993-997 (1980)

**UK & Rest of World**

184 Milton Park, Abingdon
OX14 4SE, Oxon, UK
Tel: +44 (0) 1235 828 200
Fax: +44 (0) 1235 820 482

Switzerland

Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
Tel: +41 (0) 91 604 55 22
Fax: +41 (0) 91 605 17 85

Deutschland

Bockenheimer Landsstr. 17/19
60325 Frankfurt/Main
Tel: +49 (0) 69 779099
Fax: +49 (0) 69 13376880

North America

23591 El Toro Rd, Suite #180
Lake Forest, CA 92630
Tel: + 1 800 987 0985
Fax: + 1 949 265 7703

amsbioinfo@amsbio.com

www.amsbio.com
AMS Biotechnology