

Data Sheet

JMJD2A Homogeneous Assay Kit

Catalog #AMS.50413

DESCRIPTION: The *JMJD2A Homogeneous Assay Kit* is designed to measure JMJD2A activity for screening and profiling applications. JMJD2A, also known as JHDM3A and KDM4A, is a JumonjiC (JmjC) domain containing histone lysine demethylase that exhibits demethylation activity toward H3-K₉Me³ and H3-K₃₆Me³. The *JMJD2A Homogeneous Assay Kit* comes in a convenient AlphaLISA[®] format (Scheme 1), with biotinylated histone H3 peptide substrate, primary antibody, demethylase assay buffer, and purified JMJD2A for 384 enzyme reactions. The key to the *JMJD2A Homogeneous Assay Kit* is a highly specific antibody that recognizes demethylated substrate. With this kit, only three simple steps on a microtiter plate are required for methyltransferase detection. First, a sample containing JMJD2A enzyme is incubated with the biotinylated substrate. Next, acceptor beads and primary antibody are added, then donor beads, followed by reading the Alpha-counts.

COMPONENTS:

Catalog #	Component	Amount	Storage	
AMS.50123	JMJD2A (KDM4A)	60 µg	-80°C	(Avoid freeze/thaw cycles!)
AMS.52140E	Primary antibody 5	20 µl	-80°C	
AMS.79841	Biotinylated histone H3 peptide substrate	500 rxns	-80°C	
AMS.52409	4x HDM Assay Buffer 4	3 ml	-80°C	
AMS.52031	4x Detection buffer	2 ml	-20°C	

MATERIALS REQUIRED BUT NOT SUPPLIED:


AlphaLISA[®] anti-mIgG acceptor beads, 5 mg/ml (PerkinElmer #AL105C)
 AlphaScreen[®] Streptavidin-conjugated donor beads, 5 mg/ml (PerkinElmer #6760002S)
 Optiplate -384 (PerkinElmer #6007290)
 AlphaScreen[®] microplate reader
 Adjustable micropipettor and sterile tips

APPLICATIONS: Great for studying enzyme kinetics and HTS applications.

SAFETY: This product is for research purposes only and not for human or therapeutic use. This product should be considered hazardous. Do not ingest, inhale, get in eyes, on skin, or on clothing. If so, wash thoroughly.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

Scheme 1: Our histone demethylase assays utilize highly specific antibodies that recognize demethylated products. First, a sample containing the enzyme is incubated with a biotinylated substrate. Next, acceptor beads and primary antibody are added, then donor beads, followed by reading the Alpha-counts, as shown below.

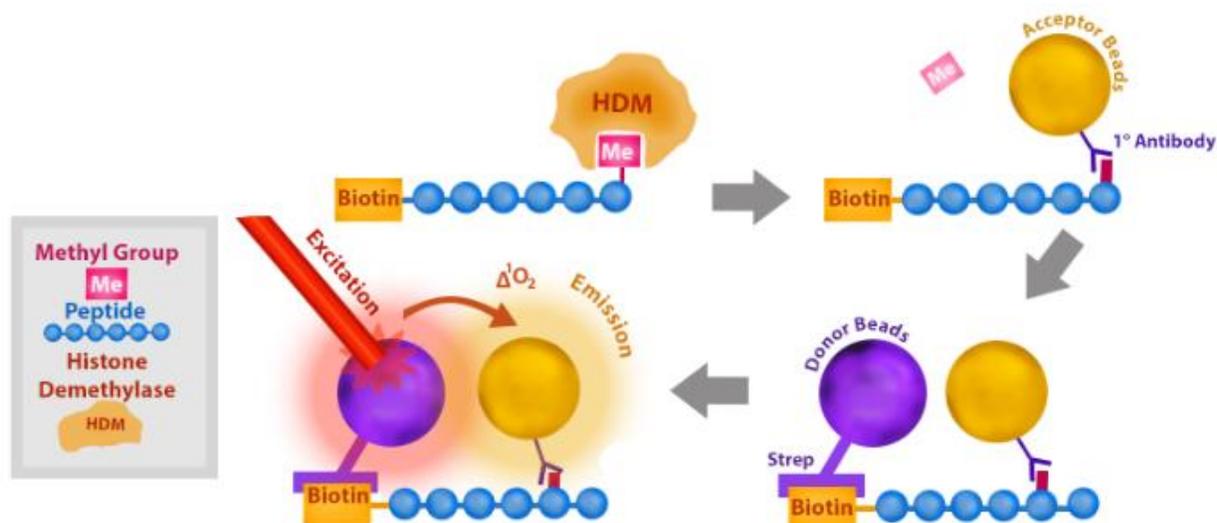

CONTRAINDICATIONS: Green and blue dyes that absorb light in the AlphaScreen signal emission range (520-620 nm), such as Trypan Blue. Avoid the use of the potent singlet oxygen quenchers such as sodium azide (NaN_3) or metal ions (Fe^{2+} , Fe^{3+} , Cu^{2+} , Zn^{2+} and Ni^{2+}). The presence of >1% RPMI 1640 culture medium leads to a signal reduction due to the presence of excess biotin and iron in this medium. MEM, which lacks these components, does not affect AlphaScreen assays.

STABILITY: At least one year from date of receipt when stored as directed.

REFERENCE(S):

1. Whetstine JR et al. *Cell* 2006; **125**: 467.
2. Fodor, B.D., *Genes & Dev.* 2006. **20**: 1557-1562.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

ASSAY PROTOCOL:

All samples and controls should be tested in duplicate. We recommend preincubating the enzyme with inhibitor, however, it is acceptable to add the substrate mixture and inhibitor followed by diluted JMJD2A without the preincubation step.

Step 1:

- 1) Re-suspend lyophilized **Biotinylated histone H3 peptide substrate** in 500 μ l of distilled water.
- 2) Prepare serial dilutions of the test inhibitors in **1x HDM Assay Buffer 4** (Scheme 2). Add 3 μ l of inhibitor solution to each well designated "Test Sample". For the wells designated "Blank" and "Positive Control" add 3 μ l of the same solution without inhibitor (typically **1x HDM Assay Buffer 4** with respective concentration of DMSO).
- 3) Thaw **JMJD2A** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full content of the tube. Aliquot **JMJD2A** enzyme into single use aliquots. Store remaining undiluted enzyme in aliquots at -80°C immediately. *Note: JMJD2A is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme.*
- 4) Dilute **JMJD2A** in **1x HDM Assay Buffer 4** at 25 ng/ μ l (100 ng/4 μ l). Keep diluted enzyme on ice until use. Discard any unused diluted enzyme after use.
- 5) Preincubate 4 μ l of diluted **JMJD2A** with 3 μ l of diluted inhibitor(s) for up to 30 minutes at room temperature, with slow shaking. For the wells designated as "Blank", add 4 μ l **1x HDM Assay Buffer 4**.
- 6) Prepare master mix: N wells \times (1.5 μ l **4x HDM Assay Buffer 4** + 1 μ l **Biotinylated substrate** + 0.5 μ l **distilled water**).
- 7) Initiate reaction by adding 3 μ l of master mix prepared as described above. Incubate at room temperature for one hour. *Note: All incubations are done with slow shaking on a rotator platform.*

Scheme 2: The serial dilution of the compounds was first performed in 100% DMSO with the highest concentration at (X) mM. Each intermediate compound dilution (in 100% DMSO) will then get directly diluted 30x fold into **1x HDM Assay Buffer 4** for 3.3x concentration (DMSO). From this intermediate step, 3 μ l of compound is added to 4 μ l of demethylase enzyme dilution is incubated for 30 minutes at room temperature. After this incubation, 3 μ l of peptide substrate is added. The final DMSO concentration is 1% for all wells.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

Reagent	Blank	Positive Control	Test Inhibitor
1x HDM Assay Buffer 4	4 µl	–	–
4x HDM Assay Buffer 4	1.5 µl	1.5 µl	1.5 µl
Biotinylated Substrate	1 µl	1 µl	1 µl
Distilled water	0.5 µl	0.5 µl	0.5 µl
Test Inhibitor/Activator	–	–	3 µl
1x HDM Assay Buffer 4 (3.3% DMSO)	3 µl	3 µl	–
JMJD2A (25 ng/µl)	–	4 µl	4 µl
Total	10 µl	10 µl	10 µl

Step 2:

Note: Protect your samples from direct exposure to light!


- 1) Dilute anti-Mouse Acceptor beads (PerkinElmer #AL105C) (1:500) and Primary antibody 5 (1:200) with 1x Detection buffer in one step. Add 10 µl of acceptor beads/antibody mixture per well. Incubate 30 min at room temperature.

Step 3:

- 1) Dilute **Streptavidin-conjugated donor beads** (PE #6760002S) 125-fold with **1x Detection buffer**. Add 10 µl of donor beads per well. Shake on a rotator platform for 30 minutes at room temperature.
- 2) Read Alpha-counts.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

Example of Assay Results:

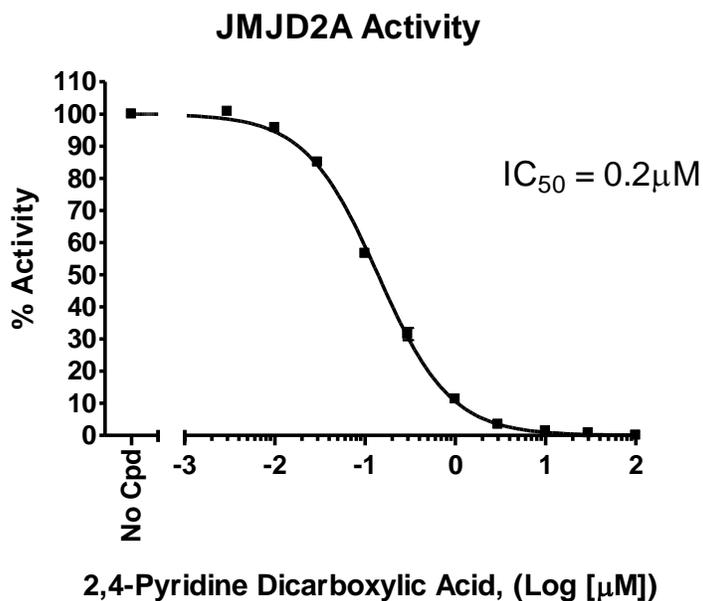

JMJD2A enzyme activity, measured using the *JMJD2A Homogeneous Assay Kit*, BPS Bioscience Cat. #AMS.50413. *Data shown is lot-specific.*

RELATED PRODUCTS:

<u>Product Name</u>	<u>Catalog #</u>	<u>Size</u>
JMJD2B Assay Kit, Homogeneous	AMS.50414-2	384 reactions
JMJD2C Assay Kit, Homogeneous	AMS.50415	384 reactions
JMJD2D Assay Kit, Homogeneous	AMS.79838	384 reactions
JMJD2E Assay Kit, Homogeneous	AMS.50417	384 reactions
JMJD2B Assay Kit, Homogeneous	AMS.50414-1	96 reactions
JMJD2C Assay Kit, Chemiluminescence	AMS.50405	96 reactions
JMJD2D Assay Kit, Chemiluminescence	AMS.50418	96 reactions
JMJD2A recombinant protein	AMS.50123	100 μg
JMJD2B recombinant protein	AMS.50111	100 μg
JMJD2C recombinant protein	AMS.50105	100 μg
JMJD2D recombinant protein	AMS.50117	100 μg
JMJD2E recombinant protein	AMS.50118	100 μg

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

AlphaScreen® and AlphaLISA® are registered trademarks of PerkinElmer, Inc.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Alpha-counts signal of positive control reaction is same as “blank” value.	JMJD2A has lost activity	Enzyme loses activity upon repeated freeze/thaw cycles. Use fresh JMJD2A, #AMS.50123. Store enzyme in single-use aliquots. Increase time of enzyme incubation. Increase enzyme concentration.
	Streptavidin Donor beads or anti-mIgG acceptor beads fail to show significant signal.	Reorder Streptavidin Donor beads or anti-mIgG acceptor beads from Perkin Elmer.
	Incorrect settings on instruments	Refer to instrument instructions for correct settings to increase sensitivity of light detection.
Alpha-counts signal is erratic or varies widely among wells	Inaccurate pipetting/technique	Run duplicates of all reactions. Use a multichannel pipettor. Use master mixes to minimize errors.

AMSBIO | www.amsbio.com | info@amsbio.com



UK & Rest of the World
184 Park Drive, Milton Park
Abingdon OX14 4SE, U.K.
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482



North America
1035 Cambridge Street,
Cambridge, MA 02141.
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218



Germany
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880



Switzerland
Via Lisano 3
(CP.683)
6900 Massagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85